

900 Touch Control



Titration was never easier


Titration at a single touch

02

Powerful as a PC, smaller than a netbook and a brilliant touchscreen – that is the 900 Touch Control. Metrohm's new benchtop control module makes titration easier, safer and faster than ever before. These are the 900 Touch Control's major functions and user benefits:

- **Touch and titrate**

Up to 14 methods can each be linked to a favorite icon on the touchscreen. You start your standard methods at a single touch. Titration cannot be easier!

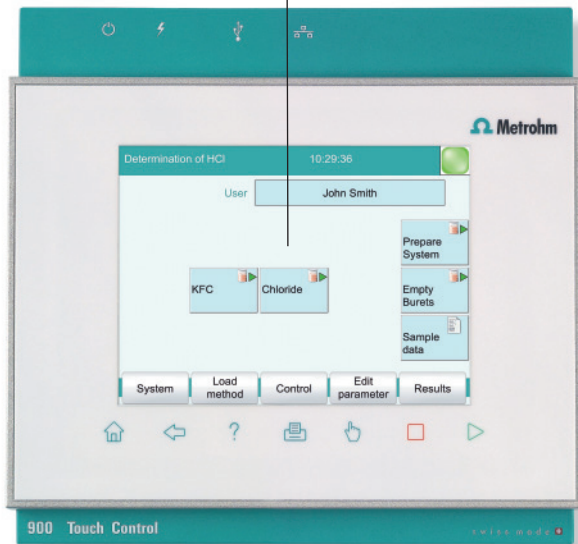

- **USB interface**

USB interface for printer, keyboard or a mouse. Simply plug in a USB stick to save your methods and determinations.



- **FDA compliance**

The 900 Touch Control complies with FDA-regulation 21 CFR, Part 11.



- **Ethernet & LIMS**

Ethernet interface for direct access to your intranet. Print reports on a network printer or save them on a network storage. Send your results directly to your LIMS or store them in the 900 Touch Control's archive.

- **Multilingual user dialog**

Ni Hao! Besides German, English, French and Spanish, the 900 Touch Control communicates also in Chinese. Further languages such as Russian and Japanese are to follow soon.

- **Automatic electrode test**

The 900 Touch Control comes with an automatic electrode test. Damaged electrodes are identified immediately; you can always rely on correct results.

- **Autostart function**

In Karl Fischer titration it can happen all too fast that sample is injected into the titration cell without the start button for the titration having been pressed first. In this case the sample is conditioned away and the determination has to be repeated. This cannot happen with the 900 Touch Control as the device starts the titration automatically once sample is added.

- **PDF Generator**

Use the 900 Touch Control's PDF generator to create forge-proof reports – without a PC!

Ordering information

900 Touch Control

2.900.0010 900 Touch Control

The 900 Touch Control is compatible with:

- 808, 809, 835, 836, 841, 842, 851, 852, 857, 888, 890, 901, 902, 904, 905, 906, 907 Titrandos
- 846 Dosing Interface
- 856 Conductivity Module
- 867 pH Module
- 814 USB Sample Processor
- 815 Robotic USB Sample Processor
- 815 Robotic Soliprep
- 855 Robotic Titrosampler
- 855 Robotic Analyzers

www.metrohm.com