



TacticID

B&W TEK

Handheld 1064nm Raman Analyzer for Rapid Identification of Suspicious Materials

The new TacticID-1064, from Delaware-based B&W Tek, is a field-ready handheld Raman analyzer designed for easy ID of a greater variety of hazardous materials, explosives, chemicals, and narcotics.

- Able to measure dark samples
- Non-destructive testing of solids, liquids, gels, and powders
- Critical component mixture analysis can identify threats in mixed samples

***Lab-Grade Results
In the Palm of Your Hand***

TacticID

Applications

- Identification of bulk explosive materials and hazardous chemicals
- On-site identification of narcotics and controlled substances
- Non-destructive forensic evidence identification
- Identification of precursors, binders, and cutting agents

Features

Ease of Use: Designed for technical and non-technical users alike, featuring simple, easy to understand menus and commands, and operation using software or hardware buttons.

Onboard Camera: Capture evidence photos to accompany scan reports.


Rugged: IP68 rated and U.S. MIL-STD-810G certified for drop testing.

Extensive Libraries: Included, with over 10,000 on-board library items - explosives*, chemicals, and narcotics with user creation and import options

Easy-to-Read Results: Color-coded Red, Yellow, and Green results with chemical safety symbols.

Bright Touch Screen: Ultra crisp HD display with touch screen capability for on-the-spot note and text inputs.

Designed for Safety: Scan delay, adjustable laser power, and hands-free operation options.

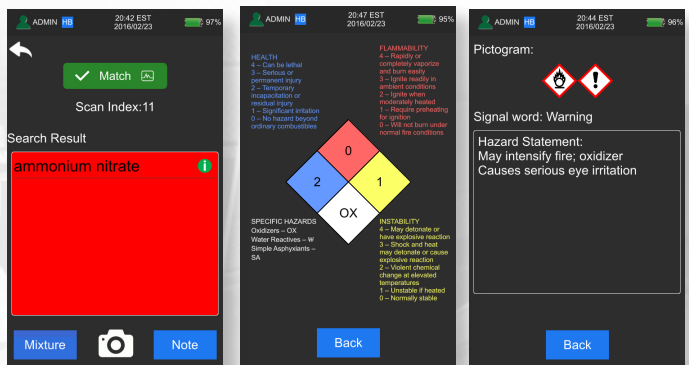

For More Information: www.tacticid.com